

MEDICAL TECHNOLOGY

We Provide Trusted & High Quality Medical Devices

"SBRG Med Tech embracing the goal by adding smiles to your life."



BRAWN GROUP OF COMPANIES is offering a better life for the people across the globe. We offer healthcare solutions at the highest level to our customer. We are accredited global healthcare company, which is empowered by one of the leading pioneer in the healthcare industry **Mr. BR Gupta,** Managing Director.





B.R. GuptaManaging Director

As a medical representative, the chasing of a big dream of serving the common man with quality medicines and anti-cancer drugs at affordable price was far from possible. What followed over the years has been the reward of sheer perseverance, passion, zeal and enthusiasm. This journey with its own ups and downs has reinforced my belief that 'Nothing is impossible', provided one religiously follows the core values- Focus, Quality, Unwavering Ethics, Transparency, Commitment and Customer Satisfaction, needless to say our strengths at **BRAWN.**



Upcoming USFDA/EU GMP

compliant plant in Halol, Gujarat



















BRAWN BIOTECH LTD.

We Are Specialized In Manufacturing & Marketing A Wide Variety of Medical Devices











Our Global Offices

SBRG Med Tech

SBRG Med Tech is committed to adding excellence in healthcare. We manufacture best medical devices. We have incorporated the latest technology and tools with the fully automated assembly systems to achieve highest quality.

Our team of experts make sure that the products offered by us fully meets the medical standards and Accuracy.

"SBRG Med Tech embracing the goal by adding smiles to your life."









MSME Certificate Certificate of Incorporation

CE-Certificate

ISO 13485:2016 Certification





Contents



■ Infusion/Transfusion Therapy

| Safety Intravenous Cannula with Injection Port | |
|---|-------|
| & with Wings | 6 |
| Intravenous Cannula with Injection Port & with Wing | s 7 |
| Intravenous Cannula | 8 |
| Extension Sets | 10-11 |
| 3 Way Stop Cock as an Accessory to Perfusion Sets | 3 12 |
| Three Way Stop Cock With Extension Line | 12 |
| Central Venous Pressure Manometer | 13 |
| Blood Administration Set / Blood Transfusion Set | 13 |
| Infusion Set for Single Use / Intravenous Infusion Se | et 14 |
| IV Flow Regulator | 14 |
| Measured Volume Set / Measured Volume IV Set | 15 |
| Scalp Vein Set | 16 |

Anaesthesia & Respiratory Care

| Suction Catheter | 17 |
|-----------------------------------|----|
| Nasal Oxygen Catheter | 18 |
| Oro-Pharyngeal Airways | 19 |
| Spinal Needle | 19 |
| Spirometer - Lungs Exerciser | 20 |
| Rebreathing Oxygen Face Mask | 20 |
| Oxygen Mask with Reservoir | 21 |
| Aerosol Therapy Mask | 21 |
| Endotracheal Tube | 22 |
| Tracheobronchial Suction Catheter | 23 |
| Ventilator Circuits | 24 |
| Breathing Circuit Bag | 25 |
| Catheter Mount | 26 |
| Heat and Moisture Exchange Filter | 26 |
| Nebulizer Mask | 27 |
| Epidural Catheter | 27 |
| Tracheostomy Tube | 28 |

Urology

| Closed-Ended Wearable Urine Collection Bag - | |
|--|----|
| with Measured Volume Meter, Sterile | 29 |
| Closed-Ended Wearable Urine Collection Bag - | |
| with Measured Volume Meter | 30 |
| Foley Catheter / Foley Balloon Catheter | 31 |
| Urine Drainage Catheters | 32 |
| Trans Uretheral Resection Set / TUR Set | 32 |
| Nelaton Catheter | 33 |

Gastroenterology

| 34 |
|----|
| 34 |
| 35 |
| 35 |
| 36 |
| 36 |
| |

Surgery and Wound Drainage

| Closed Wound Drainage System | 37 |
|--|----|
| Suction Tip and Catheter / Yankauer Suction Tube | |
| & Set with Handle | 37 |
| Trocars | 38 |
| Chest Drainage System | 39 |

Others

| Hypodermic Needle | 39 |
|---------------------------------|----|
| Umblical Artery Catheter | 40 |
| Fistula Needle | 40 |
| Surgical Dressing - I.V Fixator | 41 |
| Umbilical Occlusion Device | 41 |
| Pressure Monitoring Kit | 42 |
| | |



Safety Intravenous Cannula with Injection Port & with Wings



Cross Clip needle safety guard



Transparent flashback chamber

Ergonomically designed wings

- Extra-Sharp Back Cut bevel :- Universal back cut bevel smoothly inserts into vein with minimal puncture trauma.
- Radio Opaque catheter:- Soft catheter feels comfortable to patient and increased the indewelling time.
- Radio opaque property of catheter is good for X-rays visibility.
- Wing:- provides good stability and inclination angle of catheter to reduce stability related complications.
- Injection Port:- Allows extra passage for quick injection without interruption of main infusion line
- Needle Hub:- Ergonomic design of needle hub for better thumb support during venipuncture.
- Safety clip:- Passive safety clip prevents automatic needlestick injury.
- Flash back Chamber:- Transparent chamber to ensure quick visualization.

| | Size | Catheter O.D mm | Catheter Length (mm) | Flow Rate (ml /min) | Product Code |
|-----|--------|--------------------|-------------------------|--------------------------|-----------------|
| 14G | Orange | 2.1 | 45 | 270 | 01010614 |
| 16G | Grey | 1.7 | 45 | 200 | 01010616 |
| 17G | White | 1.5 | 45 | 125 | 01010617 |
| 18G | Green | 1.3 | 32 | 90 | 01010818 |
| 18G | Green | 1.3 | 45 | 90 | 01010618 |
| 20G | Pink | 1.1 | 32 | 60 | 01010820 |
| 22G | Blue | 0.9 | 25 | 36 | 01010922 |
| 24G | Yellow | 0.7 | 19 | 23 | 01011024 |





Intravenous Cannula with Injection Port & with Wings

Infusion/Transfusion Therapy

I.V. Cannula for continuous peripheral intravenous infusion

I.V. Cannula with port and wings

Size: 14G to 26G



| | Size | Catheter O.D mm | Catheter Length (mm) | Flow Rate (ml /min) | Product Code |
|-----|--------|--------------------|-------------------------|--------------------------|-----------------|
| 14G | Orange | 2.1 | 45 | 270 | 01010114 |
| 16G | Grey | 1.7 | 45 | 200 | 01010116 |
| 17G | White | 1.5 | 45 | 125 | 01010117 |
| 18G | Green | 1.3 | 32 | 90 | 01010318 |
| 18G | Green | 1.3 | 45 | 90 | 01010118 |
| 20G | Pink | 1.1 | 32 | 60 | 01010320 |
| 22G | Blue | 0.9 | 25 | 36 | 01010422 |
| 24G | Yellow | 0.7 | 19 | 23 | 01011524 |
| 26G | Violet | 0.6 | 19 | 15 | 01011526 |

- Special tapered & Kink resistant FEP radiopaque catheter for easy insertion with optimal flow rate
- Special flexible wings to provide easy and secure fixation
- Manufactured from tested biocompatible materials offering longer indwelling time
- Rounded grip on the needle hub for better support during insertion
- Injection port with non return silicon valve for intermittent medication




Intravenous Cannula

I.V. Cannula for continuous peripheral intravenous infusion

Intravenous Cannula without injection port & with small wings

Size: 24 G - 26G

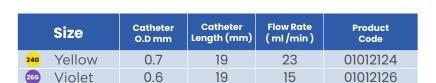

Intravenous Cannula without injection port & without wings

Size: 14G to 26G

| | Size | Catheter O.D mm | Catheter Length (mm) | Flow Rate (ml /min) | Product Code |
|-----|--------|--------------------|-------------------------|--------------------------|-----------------|
| 14G | Orange | 2.1 | 45 | 270 | 01011614 |
| 16G | Grey | 1.7 | 45 | 200 | 01011616 |
| 17G | White | 1.5 | 45 | 125 | 01011617 |
| 18G | Green | 1.3 | 32 | 90 | 01011618 |
| 18G | Green | 1.3 | 45 | 90 | 01011818 |
| 20G | Pink | 1.1 | 32 | 60 | 01011820 |
| 22G | Blue | 0.9 | 25 | 36 | 01011922 |
| 24G | Yellow | 0.7 | 19 | 23 | 01012024 |
| 26G | Violet | 0.6 | 19 | 15 | 01012026 |






Safety IV Cannula without injection port & without wings

Infusion/Transfusion Therapy

Safety I.V. Catheter without injection port and wings is designed with Passive full automatic safety feature which covers the tip of the needle completely after use, hence protecting healthcare workers against the risk of needle stick injuries and related infections.

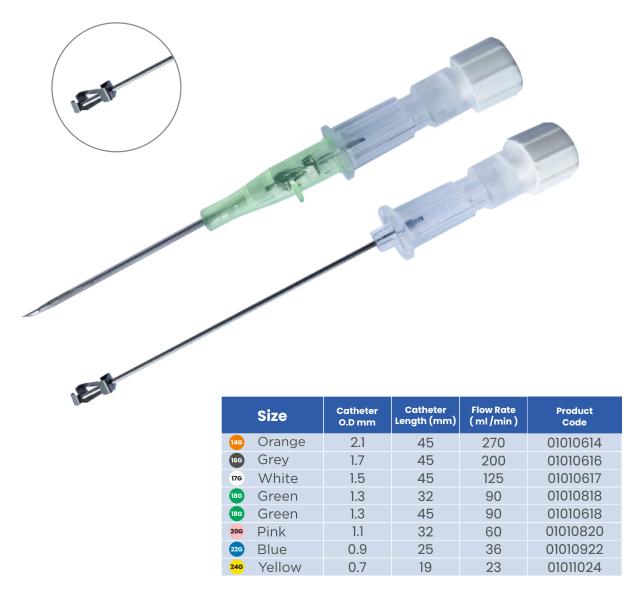


- Catheter tip is designed for a smooth cannulation experience, Radiopaque Stripes are also present for easy detection of catheter presence under X-Ray.
- Self-activation tip cover technology for protection against needle stick injuries, which cannot be bypassed.
- Universal triple beveled back cut needle, which offer lower resistance during insertion, resulting in minimum trauma.
- Gauge wise color coded body as per international standards.
- Quick Flash back confirm successful cannulation.
- DEHP / LATEX FREE.

SB NEEDLE FREE (MONO, DUO, TRIO)

Extension Sets

Extension Sets are devices that allow quick access to the infusion line, ensuring a secure connection.

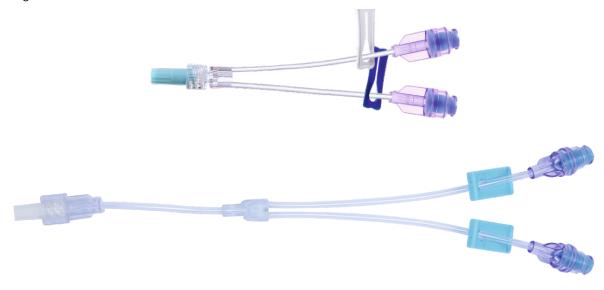

| Size | Mono | Duo | Trio | Quatro | Penta |
|--------------|----------|----------|----------|----------|----------|
| Product Code | 01070101 | 01070102 | 01070103 | 01070104 | 01070105 |

- Fast and safe access to the infusion line.
- Reduced risk of accidental puncture and infection for healthcare professionals.
- Protection of the patient's catheter from contamination.
- Reduction of the risk of device occlusion.

- Flat smooth swabbable surface for ease of disinfection.
- Available standalone and with tubing like mono, duo, trio, quatro, penta with different tube size.

| Туре | Specification |
|------------|---|
| Mono | Extension Tube small bore Mono-10 cm with needle free connector |
| Duo | Extension Tube small bore Duo-13 cm with needle free connector |
| Trio | Extension Tube small bore Trio-15 cm with needle free connector |
| Duo | Extension Tube large bore Duo-15 cm with needle free connector |
| Trio | Extension Tube large bore Trio-15 cm with needle free connector |
| Mono | Extension Tube large bore Mono-10 cm with needle free connector |
| Duo | Extension Tube small bore Duo-12 cm with needle free connector |
| Duo | Extension Tube large bore Duo-12 cm with needle free connector |
| Duo | Extension Tube small bore Duo-12 cm with needle free connector with BCV |
| Trio | Extension Tube small bore Trio-19 cm with needle free connector |
| Quatro | Extension Tube small bore Quatro-15 cm with needle free connector |
| Penta | Extension Tube small bore Penta-15 cm with needle free connector |
| Standalone | Needle free connector |



Also available as Standalone Needle Free Connector



Needle free connector to reduce the risk of Catheter Related Bloodstream Infections

EXTENSION LINES



Extension Sets

Extension Line is used for various Infusion or transfusion in order to allow for attachment with any special or general purpose medical device



High Pressure Extension Line



Low Pressure Extension Line

- Super smooth inner profile of tube ensures uniform flow rate.
- Male luer connector with cap at one end and female luer lock connector at another end.
- Available with DEHP free material.
- Transparent, Kink resistant Tubing.
- Universal 6% taper Male connection.
- Bio-compatible material
- Customized the tube length with inner & outer diameter

| Length | Product Code (EX HIGHLINE) | Product Code (EX LOWLINE) |
|--------|-------------------------------|------------------------------|
| 5 CM | 01040101 | 01040201 |
| 10 CM | 01040102 | 01040202 |
| 15 CM | 01040103 | 01040203 |
| 25 CM | 01040104 | 01040204 |
| 50 CM | 01040105 | 01040205 |
| 75 CM | 01040106 | 01040206 |
| 100 CM | 01040107 | 01040207 |
| 125 CM | 01040108 | 01040208 |
| 150 CM | 01040109 | 01040209 |
| 175 CM | 01040110 | 01040210 |
| 200 CM | 01040111 | 01040211 |
| 225 CM | 01040112 | 01040212 |
| 250 CM | 01040113 | 01040213 |
| | | |

T-WAY

Three Way Stop Cock as an Accessory to Perfusion Sets

Infusion/Transfusion Therapy

It is designed for multiple infusion of medicine, fluids and blood components to patient

Features

- Smooth Rotation of Handle align easily fluid path.
- Fully validated and can sustain pressure upto 6 bar.
- Fully transparent material allow clear visualization of fluid / medicine.
- Compatibility with other standard accessories.
- Compatible with Lipid drugs also.
- 360° rotation.



| Colour | Blue | Red | White | | |
|-----------------------------|------------------------|-----------------------|----------|--|--|
| Product Code Lipid | 01030101 | 01030102 | 01030103 | | |
| Product Code Non - Lipid | 01030201 | 01030202 | 01030203 | | |

T-WAY-EX

Three Way Stop Cock With Extension Line

It is the extension line connected one end with three way stop cock and another end with Male Luer lock, used for administration of fluids and medicines to the patients.



- Smooth Internal path without any flow turbulence.
- Fully transparent and bio-compatible material.
- Transparent, Kink resistant Tubing.
- Universal 6% taper male-female connection.

| Length | Product Code |
|--------|--------------|
| 5 CM | 01050101 |
| 10 CM | 01050102 |
| 15 CM | 01050103 |
| 25 CM | 01050104 |
| 50 CM | 01050105 |
| 75 CM | 01050106 |
| 100 CM | 01050107 |
| 125 CM | 01050108 |
| 150 CM | 01050109 |
| 175 CM | 01050110 |
| 200 CM | 01050111 |
| 225 CM | 01050112 |
| 250 CM | 01050113 |
| | |

V-PRESSCHECK



Central Venous Pressure Manometer

Infusion/Transfusion Therapy

The **Central Venous Pressure Manometer** are used for continuous or intermittent monitoring of central venous pressure during infusion

Features

- Specially designed to measure and monitor continuous or intermittent central venous pressure.
- CVP manometer is a tube with graduations marked from 4 cm to + 34 cm and attached to three way stop cock.
- Extra soft kink-resistant frosted 120cm long extension tube is fitted with female luer lock on one end and male luer lock on the other end for safe connection to central venous catheter.
- Sliding indicator is provided to record reading.
- Molded clamps with latex rubber strings are provided for easy fixation to IV Stands.
- Y-injection site/Latex rubber flash bulb is provided for extra medication.
- The manometer is manufactured from non-toxic Polypropylene & Poly Vinyl Chloride.
- The product is sterilized using EO (Ethylene Oxide) gas.

| Туре | Product Code |
|----------|--------------|
| STANDARD | 01080101 |



SB BLOOD FLOW SET

Blood Administration Set Blood Transfusion Set



- Sharp non-vented spike suitable for blood bags and all type of standard blood containers.
- Provided with a cylindrical drip chamber with 200 filter to prevent passage of any blood clots, down the path.
- Flash ball type injection port for extra medication.
- Comes with smooth 150 cm long, kink resistant tube.
- Efficient roller controller for accurate adjustment of transfusion rate.
- Sterile, individually packed.

| Туре | Product Code (B.T. Set Single Drip Chamber) | Product Code (B.T. Set Double Drip Chamber) | Product Code (B.T. Set with Y Site Single Drip Chamber) |
|------------|---|---|---|
| VENTED | 05010101 | 05010201 | 05010301 |
| NON-VENTED | 05010102 | 05010202 | 05010302 |

VSET STOP

Infusion Set for Single Use / Intravenous Infusion Set

It is used as a passage way of fluid being delivered to the patient. The transfer of the intravenous (IV) fluid is facilitated by gravitational force.

Features

- Sharp air vented spike easily puncture fluid bottle without any unwanted leakage.
- Smooth Internal path without any flow turbulance.
- Regulator, regulates the volume of flow.
- Transparent, Kink resistant Tubing.
- Universal 6% taper Male connection.
- The product is sterilized using EO (Ethylene Oxide) gas.

| Туре | Product Code |
|------------|--------------|
| Vented | 01020101 |
| Non-Vented | 01020102 |
| Auto | 01020103 |



VREG SET

IV Flow Regulator

IV Flow Regulator extension set to regulate the flow of IV fluid from an infusion set into an I.V catheter.

- Super smooth kink resistance tubing ensures uniform flow rate.
- Designed to control flow rate from 5ml/hr 250 ml/hr manually.
- Option available with "Y" connector injection site for extra medication without disturbing the main line.
- Universal Male and female luer lock makes it compatible with other devices.
- Available with DEHP free material.

| Туре | Product Code |
|------------------------|--------------|
| 'Y ' Type Connector | 01090101 |





MVF SET



Measured Volume Set / Measured Volume IV Set

Infusion/Transfusion Therapy

A burrete set is specially design to deliver a fixed volume of IV fluid at fixed rate, usually with added medication.









| Туре | Product Code |
|--------|--------------|
| 110 ml | 01060101 |
| 150 ml | 01060102 |

Features

- Specialized measured volume set made up of DEHP free material.
- Equipped with a micro drip with reduced drop size of 60 drops per ml.
- Provided with a burette type chamber of 100 ml capacity with 10 ml overflow limit.
- Burette chamber is made of bio-compatible medical grade transparent styrene polymer, suitable for infusing all types of fluids.
- Unique floating auto shut off valve acts as a floating indicator and automatically shuts off the drain path to prevent the passage of air in the line.
- No-kink device prevents kinking of tube during transportation and transfers.
- Separate plugs for extra medication and change over.
- Sterile, individually packed.

OPTIONS AVAILABLE

- Capacity : 110 ml & 150 ml
- Drop Rate : 60 drops/ml

SCALP VEIN SET

Scalp Vein Set

Infusion/Transfusion Therapy



- A short beveled siliconized needle facilitates atraumatic cannulation.
- Thin wall needle provides better flow rate per gauge.
- Butterfly shaped wings facilitate easy handling and attachment with the skin.
- Butterflies are color coded for instant identification of needle size.
- The proximal end of the set is provided with flexible female luer fitting.
- Butterfly connected to soft non-toxic, nonirritant tube which does not kink or coil.
- Sterile, individually packed in Blister pack.

| Neno Gauge | Product Code |
|------------|--------------|
| 18 | 01100118 |
| 20 | 01100120 |
| 21 | 01100121 |
| 22 | 01100122 |
| 23 | 01100123 |
| 24 | 01100124 |
| 25 | 01100125 |
| 26 | 01100126 |
| 27 | 01100127 |
| | |

SUCTION CATHETER



Suction Catheter

Anaesthesia & Respiratory Care

Suction Catheters Critically ill or injured patients who require breathing (endotracheal) tube or tracheostomy tube need occasional suctioning to remove secretions from the airway.

Thumb Control Connector





Standard Connector

Finger Tip Control Connector





Atraumatic, soft & rounded closed tip with one/two lateral eyes

| FG | Standard Connector | Thumb Control Connector | Finger Tip Control Connector | Atraumatic, Soft & rounded closed tip with one / two lateral eyes. |
|----|-----------------------|----------------------------|---------------------------------|---|
| 06 | 02010106 | 02010206 | 02010306 | 02010406 |
| 08 | 02010108 | 02010208 | 02010308 | 02010408 |
| 10 | 02010110 | 02010210 | 02010310 | 02010410 |
| 12 | 02010112 | 02010212 | 02010312 | 02010412 |
| 14 | 02010114 | 02010214 | 02010314 | 02010414 |
| 16 | 02010116 | 02010216 | 02010316 | 02010416 |
| 18 | 02010118 | 02010218 | 02010318 | 02010418 |
| 20 | 02010120 | 02010220 | 02010320 | 02010420 |
| 22 | 02010122 | 02010222 | 02010322 | 02010422 |
| 24 | 02010124 | 02010224 | 02010324 | 02010424 |

Features

- Smooth close distal end with two lateral eyes.
- Universal Female mount.
- High quality Bio-compatible material.
- Kink resistant PVC tubing with Radio-opaque line.
- Available length 50 cm.
- Size :- 5FG to 24FG.
- DEHP free material

Options Available







Light Blue







Light Green













NASAL OXYGEN CANNULA

Nasal Oxygen Catheter Twin bore nasal oxygen cannula

Nasal Oxygen Catheter is a medical device to provide supplemental oxygen therapy to people who have lower oxygen levels.



- Kink resistant medical grade tubing.
- Star lumen tube for better connection and avoid blockage.
- Efficiently administration of oxygen to the patient.
- Bio-compatible material of soft twin prong nasal tip ensure equal distribution of oxygen.

| Size | Product Code |
|------------|--------------|
| Neonatal | 04080101 |
| Paediatric | 04080102 |
| Adult | 04080103 |

ORO-PHARYNGEAL AIRWAY



Oro-Pharyngeal Airways

Anaesthesia & Respiratory Care

The **Oro-pharyngeal Airways** devices maintain un-obstucted respiration and prevent patients from biting the tongue after general anaesthesia

Features

- The oro-pharyngeal airway devices are made from the raw material of PE.
- Rounded atraumatic edges.
- Made from semi-rigid non toxic Polyethylene.
- Bite block to prevent biting of tongue and airway occlusion.
- Smooth airway path for easy cleaning.
- Color coded for easy identification of sizes
- Individually packed, EO Sterile.
- Shelf life-5 Years.
- Size: 000, 00, 0, 1, 2, 3, 4, 5



SB EASY SPINAL

Spinal Needle

Spinal Needle is used for dural puncture for both diagnostic and anesthetic purposes

- Needle hub gripping provide exceptional control during dural puncture.
- Minimize post dural puncture headache (PDPH).
- Suitable for spinal anaesthesia.
- High flow rate enables faster cerebrospinal flashback (CSF).
- Transparent needle hub helps easy visualization of CSF or blood.
- Available with customized needle length.
- Available Quincke point and Pencil point with or without introducer.



| 18 G | 19 G | 20 G | 21 G | 22 G | 23 G | 24 G | 25 G | 26 G | 27 G |
|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|
| 04090118 | 04090119 | 04090120 | 04090121 | 04090122 | 04090123 | 04090124 | 04090125 | 04090126 | 04090127 |

RESPIRATORY EXERCISER

Spirometer - Lungs Exerciser

Respiratory Exerciser is used as a comprehensive deep breathing exerciser

Features

- Restores and maintains lung capacity in post operative patients by slow, synchronized deep breathing.
- Allows visual calibration and estimation of the patients progress.
- Strengthens primary and accessory respiratory muscles and conditions them.
- Enhances endurance of both inspiratory and expiratory muscles.
- Improves respiratory parameters, Improves maximal exercise capacity.
- Improves cardio-pulmonary status of the patient, enhancing the overall fitness and well being.
- Increases circulation of hormones in the blood which increase the blood flow to the heart, brain and lungs.
- Improves oxygenation of blood, reduces fat levels by burning calories.
- Assists pulmonary rehabilitation in patients with COPD, improves the quality of their life.
- Sustained deep breathing has been shown to relieve anxiety and fight stress.



| Туре | Product Code | | |
|----------|--------------|--|--|
| STANDARD | 04100101 | | |

OXYGEN MASK

Rebreathing Oxygen Face Mask

Oxygen Mask is the device to transfer breathing oxygen gas from a storage tank to the lungs

- Smooth and feathered edge of face mask for patient's comfort and reduce irritation points.
- Super smooth kink resistance tubing ensures uniform flow rates.
- Proximal end of tube is connected with funnel shape connector for easy connection with oxygen source.
- Star lumen main tube to avoid accidental blockage.
- Tube length : 210 cm.
- Available size :- Adult, Paediatric & neonatal.



| Size | Product Code |
|------------|--------------|
| Paediatric | 04010101 |
| Adult | 04010102 |
| Neonatal | 04010103 |

OXYGEN MASK WITH RESERVOIR



Oxygen Mask With Reservoir

Anaesthesia & Respiratory Care

High Concentration oxygen mask with tubing is intended to supply oxygen or other gases to an individual for oxygen therapy, and the reservoir bag is able to avoid rebreathing.

Features

- Smooth and feathered edge of face mask for patient's comfort and reduce irritation points.
- Super smooth kink resistance tubing ensures uniform flow rate.
- Proximal end of tube is connected with funnel shape connector for easy connection with oxygen source.
- Star lumen main tube to avoid accidental blockage.
- Tube length: 210 cm.
- Available size :- Adult & Paediatric



| Size | Product Code |
|------------|--------------|
| Paediatric | 04010201 |
| Adult | 04010202 |

AEROSOL THERAPY MASK

Aerosol Therapy Mask

Aerosol Therapy Mask use oxygen, compressed air or ultrasonic power to break up solutions and suspensions into small aerosol droplets that are inhaled from the mouthpiece of the device.

| Size | Product Code |
|------------|--------------|
| Paediatric | 04020101 |
| Adult | 04020102 |

- Ideal for long term use.
- Star lumen tube for better connection.
- Uniform practical size generation (approximately 2 microns).
- Very effective for medication delivery.
- Light weight and easy to handle.




ENDOTRACHEAL TUBES

Endotracheal Tube

It is used to provide oxygen and inhaled gases to the lungs and protects the lungs from contamination, such as gastric contents or blood.



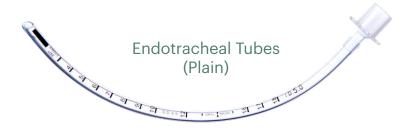

- High Volume, low pressure cuff, maintains seal throughout ventilator cycle.
- X-ray detective line available.
- Thermosensitive PVC allows the tube to soften at body temperature to confirm airway anatomy
- Clear PVC tube with Visible marking for easy observation.
- Murphy eye smoothly formed to allow ventilation in the event of obstruction to the end of tube during intubation.
- Universal 15mm Connector at proximal end.
- Available with DEHP free Material.
- Sterile/Disposable / Individually packed.

| Size/MM | PLAIN ET TUBE PRODUCT CODE | CUFFED ET TUBE PRODUCT CODE |
|---------|-------------------------------|--------------------------------|
| 2.0 | 04030101 | 04030201 |
| 2.5 | 04030102 | 04030202 |
| 3.0 | 04030103 | 04030203 |
| 2.5 | 04030104 | 04030204 |
| 4.0 | 04030105 | 04030205 |
| 4.5 | 04030106 | 04030206 |
| 5.0 | 04030107 | 04030207 |
| 5.5 | 04030108 | 04030208 |
| 6.0 | 04030109 | 04030209 |
| 7.5 | 04030110 | 04030210 |
| 8.0 | 04030111 | 04030211 |
| 8.5 | 04030112 | 04030212 |
| 9.0 | 04030113 | 04030213 |
| 9.5 | 04030114 | 04030214 |
| 10.0 | 04030115 | 04030215 |
| | | |

SB BREATHING CATHETER



Tracheobronchial Suction Catheter

Anaesthesia & Respiratory Care

Tracheobronchial Suction Catheter is a device which is indicated for use in new-borns, infants or small children to clear the secretion from oral cavity.



- Specially designed for aspiration of secretions from oropharynx in newly born babies.
- Clear transparent container permits immediate visual examination of the aspirate.
- Also suitable for obtaining mucus specimen for micro biological examinations.
- Spare plug cap is provided to seal the container for safe transportation of specimen to laboratory or aseptic disposal of container.
- Sterile, individually packed in a peelable medical grade paper pack / poly bag

| FG | Product Code |
|----|--------------|
| 08 | 04060108 |
| 10 | 04060110 |
| 12 | 04060112 |
| 14 | 04060114 |


SB VENT CIRCUIT


Ventilator Circuit

Ventilator circuit is the gas pathway connecting patients lung to the mechanical ventilator.



Ventilator Circuits with one water trap

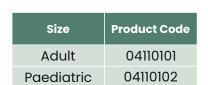

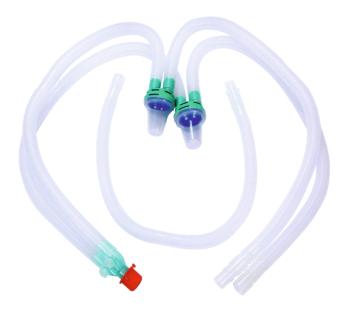

Ventilator Circuits with two water trap

| Size | Product Code |
|------------|--------------|
| Adult | 04110201 |
| Paediatric | 04110202 |

- Re- sealable water-trap provides a perfect sealing while removal of excess collection of water from the water-traps giving a continuous flow.
- Made of high quality medical grade plastic with 22M-22F straight connector.
- Transparent for easy monitoring of water build-up or blockages.
- Corrugated, Flexible and Light weight tubing 1.6 meter with single and double re-sealable water trap, 'Y Piece' with monitoring port.
- Individually packed in poly pack, Available sizes : Adult and Pediatric.

SB BREETH CIRCUIT



Breathing Circuit Bag

Anaesthesia & Respiratory Care

Bain circuits is a co-axial modification of the basic T-piece system, developed to facilitate scavenging of waste anesthetic gases.



- Compact and inexpensive, Low resistance to breathing.
- Facilitates scavenging of waste gases
- Adult :- It is based on Mapleson's D system, Provided with 22 mm corrugated tube, expiratory valve and 2 litre antistatic bag.
- Pediatric:- It is based on Mapleson's F system, provided with 15 & 10 mm corrugated tube with a 0.5 liter antistatic rebreating bag and a bag bleed valve.
- Sterile/ Disposable / Individually packed.

| Туре | Product Code |
|----------|--------------|
| STANDARD | 04120101 |

CATHETER MOUNT


Catheter Mount

Catheter Mount connects patient ET to breathing circuit. This minimizes the jerk of breathing to ET or TT reducing trauma to trachea.

Features

- The double swivel makes no dragging on breathing circuit & ET or TT.
- Light weight double swivel rotates 360° on both axis, minimum compliance, Low dead.
- Compatible with all type of breathing and ventilator circuits.
- Provided with 22 mm standard female connector on both the ends with inner diameter of 15 mm.
- Collapsible corrugated tubing being inert to all anesthetic gases and all reagents.
- An elbow port is provided for easy suctioning without causing disturbance to breathing circuits.
- Sterile/ Disposable / Individually packed.

| Туре | Product Code |
|----------|--------------|
| STANDARD | 04070101 |

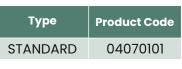

HME FILTER

Heat and Moisture Exchange Filter

HME Filters used to maintain passively heat and moisturizer levels of air prevent any bacteria or viruses from entering the patient's airways.

- Provides effective protection against various types of particles including bacteria, viruses, and moisture droplets.
- Precisely engineered and easy to use device made up of high strength plastic material offering reliable performance.
- Provided with high efficiency electrostatic filters to address concerns associated with ventilated patients.
- The filters help protect the patient, the equipment, and the breathing circuit from contamination.
- Enhances particle capture efficiency and give a low and stable breathing resistance even when wet.
- A conveniently located port or connection permits gas sampling at the top of the housing,
- H.M.E.F's are designed according to "ISO" standards to help facilitate connection to different ventilators and anesthesia systems.





| Size | Product Code |
|------------|--------------|
| Paediatric | 04040101 |
| Adult | 04040102 |

NEBULIZER MASK



Nebulizer Mask

Anaesthesia & Respiratory Care




Features

- Soft, clear Aerosol Mask with anatomical form is ideal for long term use.
- Gently rolled, feathered edges with integrated nose bridge for extra comfort.
- Swivel connector for convenience of nebulization in horizontal and vertical position.
- Larger surface area provided by unique convex cone design, ensures maximum capillary action, and eliminates medicinal wastage.

| Size (210 cm Tube) | Product Code |
|-----------------------|--------------|
| Child | 04020101 |
| Adult | 04020102 |

SB EPI ANESTHESIA

Epidural Catheter

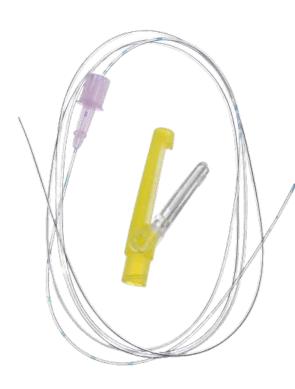

- Kink and stretch resistant.
- Easy to thread and remove.
- Latex-free.
- Easy to anchor with the Catheter Connector
- The catheter has marking every centimeter from 5 to 20 cm and black mark at the distal end to ensure that entire catheter has been withdrawn.
- A catheter guide which connects on the Tuohy needle facilitates the insertion of the catheter in the epidural space.
- Yellow label for epidural access identification.

| Gauge | Product Code |
|-------|--------------|
| 19 | 04140119 |
| 20 | 04140120 |

SB FLEXI TRACHEAL TUBE

Tracheostomy Tube

Anaesthesia & Respiratory Care



- Made of thermosensitive material which softens at body temperature to accommodate individual patient anatomy.
- Soft flexible flange helps easy fixation.
- 15 mm connector can be rotated in 360-degree direction & facilitates the safe connection to breathing system from any direction.
- Radio opaque blue line for X-ray visualization.
- Soft low-pressure cuff dynamically adjusts the pressure at the point of contact during the breathing cycle. Thus minimize the risk of tracheal trauma & other complications.
- Pilot balloon with non-return valve facilitates the cuff inflation & deflation.
- Supplied with cotton neck tape.
- Sterile, individually Blister packed.

URI COLLECTO



Closed-Ended Wearable Urine Collection Bag with Measured Volume Meter, Sterile

Urology

A **Urine Collection** device allows an individual to empty their bladder into a container hygienically and without spilling urine.



| Size | Product Code |
|------------|--------------|
| Paediatric | 03010101 |
| Adult | 03010102 |



- Super smooth kink resistance tube provided with universal tapered connector.
- Specially designed hook or hanger facilities for carrying, handling and holding the tube in upright position.
- Efficient non-return valve prevents the back flow.
- Bag graduated in ml to indicate the quantity of urine collected.
- Transparent sheeting allows visual inspection.

URI COLLECTO MV

Closed-Ended Wearable Urine Collection Bag - with Measured Volume Meter

Urine collection bag graduated in ml with measured volume meter for accurate measurement of urine

| Product Code |
|--------------|
| 03010201 |
| 03010202 |
| |





Features

- Super smooth kink resistance tube provided with universal tapered connector.
- Specially designed hook or hanger facilities for carrying, handling and holding the tube in upright position.
- Designed for accurate measurement of hourly urine output during examination of the patient.
- Urine bag attached directly with measured volume meter ensures auto over flow and convenience to empty the measured volume meter in the bag.
- Efficient non-return valve prevents the back flow.
- Bag graduated in ml to indicate the quantity of urine collected.
- Transparent sheeting allows visual inspection.
- With "T" type bottom outlet.

OPTIONS AVAILABLE

- With or without pinch clamp.
- Chamber size: 250 ml, 500 ml.
- Tube Length :- 100 cm, 120 cm, 150 cm.
- Bag size :- 2000 ml.
- With or without sampling port.

P-DRAIN BALLOON CATH



Foley Catheter / Foley Balloon Catheter

Urology

Foley balloon is a flexible tube that a clinician passes through the urethra and into the bladder to drain urine



| FG | Two way Foley balloon catheter (Latex) | Three way Foley balloon catheter (Latex) | Two way Foley balloon catheter (Silicon) | Three way Foley balloon catheter (Silicon) |
|----|--|--|--|--|
| 06 | 03020106 | 03020206 | 03020306 | 03020406 |
| 08 | 03020108 | 03020208 | 03020308 | 03020408 |
| 10 | 03020110 | 03020210 | 03020310 | 03020410 |
| 12 | 03020112 | 03020212 | 03020312 | 03020412 |
| 14 | 03020114 | 03020214 | 03020314 | 03020414 |
| 16 | 03020116 | 03020216 | 03020316 | 03020416 |
| 18 | 03020118 | 03020218 | 03020318 | 03020418 |
| 20 | 03020120 | 03020220 | 03020320 | 03020420 |
| 22 | 03020122 | 03020222 | 03020322 | 03020422 |
| 24 | 03020124 | 03020224 | 03020324 | 03020424 |
| 26 | 03020126 | 03020226 | 03020326 | 03020426 |

- Non-toxic, bio-compatible and extra smooth for maximum patient comfort.
- Symmetrical balloon, rounded sealed tip and hard non-return.
- X-ray opaque line allows for confirmation of intubated tube using X-ray.
- Soft and uniformly inflated balloon makes the tube sit well against the bladder.
- Smooth round shaft can minimize trauma during insertion and withdrawal.
- Comes as 2-way or 3-way design. 3-way design helps in irrigation of bladder and sometimes for therapeutic instillation.
- Colour coded for easy identification.

URINE DRAINAGE CATHETER

Urine Drainage Catheters

They are intended to be inserted intermittently through the urethra and up to the bladder enabling urine to drain. They are designed for short-term, intermittent transient use only.





Soft Tip



Features

- Super smooth kink resistance tube provided with universal tapered connector.
- Specially designed hook or hanger facilities for carrying, handling and holding the tube in upright position.
- Efficient non-return valve prevents the back flow.
- Bag graduated in ml to indicate the quantity of urine collected.
- Transparent sheeting allows visual inspection.


| FG | Product Code |
|----|--------------|
| 06 | 03020606 |
| 80 | 03020806 |
| 10 | 03021006 |
| 12 | 03021206 |
| 14 | 03021406 |
| 16 | 03021606 |
| 18 | 03021806 |
| 20 | 03022006 |
| 22 | 03022206 |
| 24 | 03022406 |

T-URETHSET

Trans Uretheral Resection Set TUR Set

TUR Set is used for endoscopic irrigation during the Trans Urethral Resection of the prostate gland.

- Manufactured from clinical-grade non-toxic, transparent medical PVC sheet.
- Y shaped connector facilitates alternative changes solution.
- Thumb operated clamps smooth change of bottle.
- Proximal end fitted with flexible latex tubing for easy connection to the endoscope.
- ETO Sterile, individually packed.

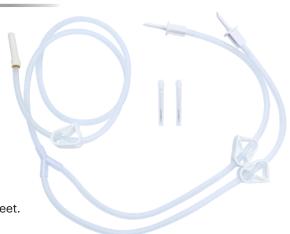

NELATON CATHETER



Nelaton Catheter Urology



| FG | (20cm) | (40 cm) |
|----|----------|----------|
| 06 | 03020506 | 03020606 |
| 80 | 03020508 | 03020608 |
| 10 | 03020510 | 03020610 |
| 12 | 03020512 | 03020612 |
| 14 | 03020514 | 03020614 |
| 16 | 03020516 | 03020616 |
| 18 | 03020518 | 03020618 |
| 20 | 03020520 | 03020620 |
| 22 | 03020522 | 03020622 |
| 24 | 03020524 | 03020624 |

- Nelaton catheters are used for short term bladder catheterization.
- Manufactured from non-toxic, non-irritant medical grade PVC compatible with catheter lubricants.
- Frozen surface tubing for super smooth intubation.
- The presence of radio opaque line throughout the catheter tube.
- The distal end is coned for non-traumatic introduction with two lateral eyes for efficient drainage.
- The proximal end is fitted with universal funnel shaped connector for extension.
- Funnel shape connectors are color coded for instant identification of sizes.
- Straight packed in peelable soft blister pack.
- Sterile, ready for use.

RYLE'S TUBE

Ryle's Tube

Ryle's tube is a narrow-bore tube passed into the stomach via the nose. It is used for short- or medium-term nutritional support, and also for aspiration of stomach contents - eg, for decompression of intestinal obstruction.

Features

- Smooth close distal end with two lateral eyes.
- Universal Female mount.
- High quality Bio-compatible material.
- Kink resistant PVC tubing with Radio-opaque line.
- Available length 40, 50 & 100 cm.
- DEHP free material
- Size :- 4 FG to 10 FG.





| FG | Product Code |
|----|--------------|
| 08 | 02020108 |
| 10 | 02020110 |
| 12 | 02020112 |
| 14 | 02020114 |
| 16 | 02020116 |
| 18 | 02020118 |

INFANT FEEDING TUBE

Feeding Tube

It is used to feed a person who is unable to eat or drink. A person may need a feeding tube due to difficulty swallowing, an eating disorder, or other feeding issues.

Features

- Smooth close distal end with two lateral eyes.
- Universal Female mount.
- High quality Bio-compatible material.
- Kink resistant PVC tubing with Radio-opaque line.
- Available length 40, 50 & 100 cm.
- Size :- 4 FG to 10 FG.
- DEHP free material.



Female Mount

| Size | Product code Feeding tube (40 cm) | Product code Feeding tube (50 cm) | Product code Feeding tube (100 cm) |
|------|---|---|--|
| 04 | 02030104 | 02030204 | 02030304 |
| 06 | 02030106 | 02030206 | 02030306 |
| 80 | 02030108 | 02030208 | 02030308 |
| 10 | 02030110 | 02030210 | 02030310 |

FEEDING BAG



Feeding Bag

Gastroenterology

It is used to feed a person who is unable to eat or drink. A person may need a feeding tube due to difficulty swallowing, an eating disorder, or other feeding issues.

Features

- Specially designed 1.2 Lit. capacity bag for feeding purpose.
- An ice pouch is provided with the bag for preserving the food for longer period.
- A transfer set is directly connected to the bag for easy feeding.
- Roller controller is provided to control the flow at which food is going to the patient.
- A hanger is provided with bag for easy handling.
- Sterile, individually packed in peelable pouch.





UMBILICAL CATHETER

Umblical Artery Catheter

An umbilical artery catheter (UAC) allows blood to be taken from an infant at different times, without repeated needle sticks. It can also be used to continuously monitor a baby's blood pressure.





- A traumatic, soft rounded open distal tip to reduce vascular trauma.
- Tube with radio-opaque line, marked at every cm from 5 to 25 cm from the open distal tip for accurate placement
- Can be used with venous or arterial routes for infusion, transfusion, administration of medication, blood sampling, pressure measurement etc.
- Available with female flexible mount and luer lock.
- Color coded funnel end connector for easy identification of size.
- Length: 40 cm, Sizes: 4, 5, 6 & 8 FG



| Size/Colour | Product Code |
|-------------|--------------|
| 4G Red | 02040104 |
| 5G Grey | 02040105 |
| 6G Green | 02040106 |
| 8G Blue | 02040108 |

STOMACH EVACULATOR

Stomach Evaculator

Gastroenterology



| FG | Product Code |
|----|--------------|
| 06 | 02050106 |
| 08 | 02050108 |
| 10 | 02050110 |
| 12 | 02050112 |
| 14 | 02050114 |
| 16 | 02050116 |
| 18 | 02050118 |
| 20 | 02050120 |
| 22 | 02050122 |
| | |

Features

- Suitable for nasogastric introduction, for nutritional purposes or aspiration of intestinal secretions.
- The distal end is coned with corrosion resistant stainless-steel balls sealed into the tube.
- Manufactured from medical grade Polymer.
- Radio opaque line is provided throughout the length.
- Marked at 50, 60, 70 cm: Sterile, individually packed.

LEVINE TUBE

Levine Tube



- Features
- Material: non-toxic medical PVC.
- Smooth and transparent tube (or frosted tube).
- Rounded opened distal end.
- Marked at intervals from tip for accurate placement.
- Radio opaque line throughout the length for X-ray

| FG | Product Code |
|----|--------------|
| 06 | 02060106 |
| 08 | 02060108 |
| 10 | 02060110 |
| 12 | 02060112 |
| 14 | 02060114 |
| 16 | 02060116 |
| 18 | 02060118 |
| 20 | 02060120 |
| 22 | 02060122 |
| | |

SB CLOSED CARE DRAIN



Closed Wound Drainage System

Surgery and Wound Drainage

It is used to reduce postoperative edema and hematoma formation, decrease the possibility of infection, and minimize the chance of external contamination of the surgical site.

Features

- Sizes for accessories: 8, 10, 12, 14, 16 & 18 FG
- Bellow unit with connector
- Connecting tube with clamp and "Y" connector
- Curved trocar with catheter
- Spare perforated Redon drain catheter.
- Bellow size :- 600 ml & 800 ml

| Size (FG) | Product code Closed Wound Suction Units (600 cm) | Product code Closed Wound Suction Units (800 cm) |
|-----------|--|--|
| 08 | 06010108 | 06010208 |
| 10 | 06010110 | 06010210 |
| 12 | 06010112 | 06010212 |
| 14 | 06010114 | 06010214 |
| 16 | 06010116 | 06010216 |
| 18 | 06010118 | 06010218 |
| | | |

SUCTION TIP AND CATHETER

Suction Tip and Catheter / Yankauer Suction Tube & Set with Handle

Suction Tip and Catheter is a suction set used for the removal of blood and fluids during surgery

Type

Plain Tip

Crown Tip





Standard Tip

Product Code

06020101

06020102


- Tube length is 250 cm.
- Soft, flexible & kink resistant tube
- Ribs on the outer surface for better grasping.
- Siliconized sleeve on the handle for smooth operation.
- Suitable for intra-surgical removal of secretion and body fluids.
- Thumb control Yankauer handle available in standard or crown tip.



THORACIC DRAINAGE CATHETER

Trocars

It is used to remove air in the case of pneumothorax or fluid such as in the case of pleural effusion, blood, chyle, or pus when empyema occurs from the intrathoracic space

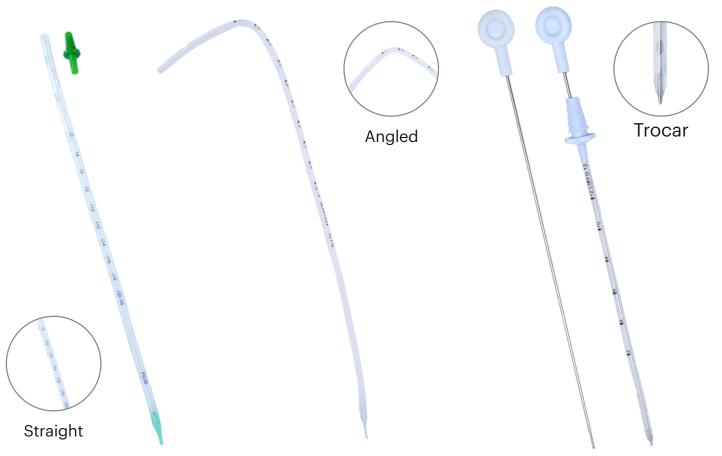

- Super Smooth Kink Reistance tubing ensures uniform flowrate.
- Catheter Without Trocar:- Atraumatic, Soft Rounded, Open Distal end with six lateral eyes for non-traumatic insertion.
- Catheter with Trocar:- Atraumatic, soft rounded, open distal end with two lateral eyes for non-traumatic insertion.
- Cross side eyes to prevent tissue aspiration.
- Proximal end is fitted with tapered connector for easy connection to drainage bottle.
- X-ray opaque line provided throughout the length of catheter.
- Marking on catheter is provided at from the distal end to help ascertain placement.
- Available with DEHP free material.

| Size | Thoracic Drainage Catheter Straight Product Code | Thoracic Drainage Catheter Angled Product Code | Thoracic Drainage Catheter Trocar Product Code |
|------|--|--|--|
| 08 | 06030108 | 06030208 | 06030308 |
| 10 | 06030110 | 06030210 | 06030310 |
| 12 | 06030112 | 06030212 | 06030312 |
| 14 | 06030114 | 06030214 | 06030314 |
| 16 | 06030116 | 06030216 | 06030316 |
| 18 | 06030118 | 06030218 | 06030318 |
| 20 | 06030120 | 06030220 | 06030320 |
| 22 | 06030122 | 06030222 | 06030322 |
| 24 | 06030124 | 06030224 | 06030324 |
| 26 | 06030126 | 06030226 | 06030326 |
| 28 | 06030128 | 06030228 | 06030328 |
| 30 | 06030130 | 06030230 | 06030330 |
| 32 | 06030132 | 06030232 | 06030332 |
| 34 | 06030134 | 06030234 | 06030334 |
| 36 | 06030136 | 06030236 | 06030336 |
| 38 | 06030138 | 06030238 | 06030338 |
| 40 | 06030140 | 06030240 | 06030340 |
| | | | |

SB DRAIN KIT



Chest Drainage System / Under Water Seal Drainage System

Surgery and Wound Drainage

Underwater seal **Chest Drainage** is used to restore proper air pressure to the lungs, reinflate collapsed lungs as well as remove blood and other fluids.

Features

- Moderate to heavy drainage.
- Accurate, graduated marked on the bottle to determine the amount of collection.
- Clearly marked initial level ensures effective under water seal.
- Complete positive and negative pressure controls.
- Proprietary unit for pneumonectomy.
- Color-coded, needleless Quick-disconnect tubing for patient blood sampling.
- Generous collection chambers holding up to 500, 1200 & 2000 cc of fluid.
- Double chamber and single chamber compact unit suitable for all cardiac and thoracic procedures.
- Range appropriate for controlling massive, moderate, and small drainage systems.
- There is a suction port available for connecting to the suction device.
- Large bore kink resistant tubing allows for unhindered flow.
- Hangers made of metal for simple connection to the patient's bedside.
- Sterile, individually packed.



HYPODERMIC NEEDLE

Hypodermic Needle

Others



- A proprietary grinding process results in an exceptionally sharp lancet point able to penetrate tissue smoothly and with minimal trauma.
- Hub is color-coded for accurate, quick identification in or out of the package.
- Hub translucency confirms flashback.
- Paper-to-plastic packaging is convenient and environmentally friendly.
- Available in gauges 18-27 and 30 sizes.

| Gauge | Product Code |
|-------|--------------|
| 18 | 08050118 |
| 19 | 08050119 |
| 20 | 08050120 |
| 21 | 08050121 |
| 22 | 08050122 |
| 23 | 08050123 |
| 25 | 08050125 |
| 27 | 08050127 |
| 30 | 08050130 |

UMBILICAL ARTERY CATHETER

Umblical Artery Catheter

Features

- Designed for intermittent or continuous access to the umbilical vessels of newborn or pre-mature infants.
- Suitable for monitoring blood pressure, withdrawal of blood for sampling without repeated needle pricks, intra venous
- Infusion and exchange transfusion etc.
- Smooth round tip helps easy cannulation and silk finish surface further facilitates the passage of catheter in the vein / artery.
- Open distal end without lateral eyes, eliminates the chance of clot formation in the blind spaces.
- Radio opaque line is provided for X-ray visualization.
- Clear scale markings at every cm helps to ascertain depth of insertion.
- Proximal end is fitted with female luer mount for convenient connection to I.V therapy equipment.
- Individually packed in peelable Pack.
- Sterile, ready for use.


| FR | Product Code | |
|-----|--------------|--|
| 3.5 | 08030135 | |
| 4.0 | 08030140 | |
| 5.0 | 08030150 | |
| 6.0 | 08030160 | |
| 8.0 | 08030180 | |



FISTULA NEEDLE

Fistula Needle

- Back eye for uninterrupted blood flow through the needle.
- Ratable wings for firm grip, easy turning, and secure fixation.
- Smooth Geometry for laminar blood flow.
- Modified hub design with black mark for positioning the bevel accurately.
- Kink resistant tubing with smooth occlusion clamp.
- Color coded wings for identification for needles size.
- Option for single/Twin Fistula needles.
- Disposable, sterile, non-pyrogenic.

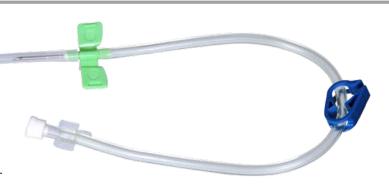

| Gauge | Product Code | |
|-------|--------------|--|
| 15 | 08040115 | |
| 16 | 08040116 | |
| 17 | 08040117 | |

SB CANN FIX



Surgical Dressing - I.V Fixator

Others

Cannula Fixator with extra strap for additional reinforcement around the cannula. It also prevents accidental movement of the catheter and prevent infection.

Features

- Used for fixation of I.V. Cannula / Catheter.
- Made of thin Polyurethane film prevent contaminating the infusion site.
- Breathable and water proof, with reinforced notches that reduce mechanical stress from heavy catheters.
- Thin film ensures easy fixation of devices with contoured surfaces.
- Available in two variants Transparent and Non-woven.
- Sterile / Disposable / Individually packed





| Туре | Product Code | |
|----------|--------------|--|
| STANDARD | 07010101 | |

UC CLAMP

Umbilical Occlusion Device

It is used for binding of the umbilical cord by nipper to interrupt blood flow from placenta to foetus.









- Designed finger grip provision for easy to use.
- High quality Bio-compatible material.
- Serration along the length prevent slipping of the umbilical cord.
- Manufactured from non toxic medical grade polymer.
- Easy Re-open facilities.

| Туре | Product Code | |
|----------|--------------|--|
| STANDARD | 07020101 | |

PRESSURE MONITORING KIT

Pressure Monitoring Kit



| ТҮРЕ | SINGLE | DOUBLE | TRIPLE |
|--------------|----------|----------|----------|
| Product Code | 08060101 | 08060102 | 08060103 |

- Pressure Monitoring Kits are a combination of medical grade PVC Tubing, Drip Chambers, Transducers, Flush device & Stop cocks, sampling devices, connectors, roller clamps & caps assembled as a functional unit to monitor the invasive blood pressure.
- This functional unit is sterilized by Ethylene Oxide, as a "single unit sterile pack".
- The Sensor is a hydraulic communication with the blood stream to detect the patient's pulse. The flush device allows continuous flow of 3 ml/hr. When addressed when operated it allows high flow, that is required for flushing blood remains from the monitoring line. The stopcock is used for Zeroing the pressure monitor before operation. It also assists in Blood samplings procedure.
- Pressure transducers are placed on the base plate which is usually mounted on IV stand. Transducers should be positioned at the level of the patient's heart to get accurate reading.
- The whole apparatus is primed with saline solution & then connected to the site of monitoring (Radial Artery in the wrist via a catheter or cannula).
- Saline is connected into the pressure bag & inflated up to 300 mm hg to prevent backflow of blood from monitoring site.





SBRG Med Tech

embracing the goal by adding smiles to your life.





MEDICAL TECHNOLOGY- PRODUCT CATALOGUE

Plot No.-352, Pace City-2, Sector 37, Gurugram - 122001, Haryana, INDIA Email: info@sbrgmedtech.in | regulatory@sbrgmedtech.in

FARIDABAD OFFICE: 13, NIT, Industrial Area, Faridabad-121001, Haryana, India

CORPORATE OFFICE: Plot No 30 Sector 33 Near Hero Honda Factory Gurgaon- 122001 (Haryana) India DELHI HEAD OFFICE: A-36, First Floor, Guru Nanak Pura, Laxmi Nagar, East Delhi, Shahdara, Delhi-110092